



**Europäisches  
Patentamt**

**European  
Patent Office**

**Office européen  
des brevets**

RECEIVED

14 APR 2004

WIPO

PCT

**Bescheinigung**

**Certificate**

**Attestation**

Die angehefteten Unterla-  
gen stimmen mit der  
ursprünglich eingereichten  
Fassung der auf dem näch-  
sten Blatt bezeichneten  
europäischen Patentanmel-  
dung überein.

The attached documents  
are exact copies of the  
European patent application  
described on the following  
page, as originally filed.

Les documents fixés à  
cette attestation sont  
conformes à la version  
initialement déposée de  
la demande de brevet  
européen spécifiée à la  
page suivante.

**Patentanmeldung Nr. Patent application No. Demande de brevet n°**

03005186.6

Best Available Copy

Der Präsident des Europäischen Patentamts;  
Im Auftrag

For the President of the European Patent Off

Le Président de l'Office européen des brevet  
p.o.

**R C van Dijk**

**PRIORITY  
DOCUMENT**

SUBMITTED OR TRANSMITTED IN  
COMPLIANCE WITH RULE 17.1(a) OR (b)



Anmeldung Nr:  
Application no.: 03005186.6  
Demande no:

Anmeldetag:  
Date of filing: 07.03.03  
Date de dépôt:

Anmelder/Applicant(s)/Demandeur(s):

ISTITUTO NAZIONALE PER LO STUDIO  
E LA CURA DEI TUMORI  
Via Venezian, 1  
I-20133 Milano  
ITALIE

Bezeichnung der Erfindung/Title of the invention/Titre de l'invention:  
(Falls die Bezeichnung der Erfindung nicht angegeben ist, siehe Beschreibung.  
If no title is shown please refer to the description.  
Si aucun titre n'est indiqué se referer à la description.)

Anaplastic limphoma kinase assay, reagents and compositions thereof

In Anspruch genommene Priorität(en) / Priority(ies) claimed /Priorité(s)  
revendiquée(s)  
Staat/Tag/Aktenzeichen/State/Date/File no./Pays/Date/Numéro de dépôt:

Internationale Patentklassifikation/International Patent Classification/  
Classification internationale des brevets:

G01N33/48

Am Anmeldetag benannte Vertragstaaten/Contracting states designated at date of  
filing/Etats contractants désignées lors du dépôt:

AT BE BG CH CY CZ DE DK EE ES FI FR GB GR HU IE IT LU MC NL  
PT SE SI SK TR LI RO

## ANAPLASTIC LYMPHOMA KINASE ASSAY, REAGENTS AND COMPOSITIONS THEREOF

The present invention provides an assay for measuring the kinase activity of the Anaplastic lymphoma kinase (ALK). More specifically the invention is directed to methods, reagents and compositions for detecting the tyrosine-phosphorylating activity of ALK. The assay is particularly useful for  
5 in vitro screening and identification of potential ALK inhibitors.

### **BACKGROUND OF THE INVENTION**

Anaplastic lymphoma kinase (ALK) is a receptor tyrosine kinase, believed to play an important role in the development and function of the nervous system. ALK is normally expressed in the central nervous system,  
10 with peak expression during the neonatal period. However, due to chromosomal translocations, ALK is also aberrantly expressed and activated in some cancers in the form of oncogenic fusion proteins. ALK fusion proteins are responsible for approximately 5 - 10% of all non-Hodgkin's lymphomas. The annual incidence of ALK positive lymphomas is about 100 000  
15 worldwide, with 2000 - 3000 new cases occurring in EU countries. ALK is an excellent candidate for therapeutic intervention, as it plays an essential role in oncogenicity and its normal expression is mostly restricted to the central nervous system. Hence, a specific ALK inhibitor could be an efficient treatment for ALK positive lymphomas with few associated clinical side  
20 effects. For the identification of potential ALK inhibitors it is highly desirable to develop a suitable assay for assessing directly the inhibitory effect of compounds on enzyme activity.

### **DESCRIPTION OF THE INVENTION**

According to the invention, an *in vitro* kinase assay specific for ALK  
25 and useful for screening compounds that modulate ALK activity is provided.

This assay is based on the use of a peptide substrate which is readily phosphorylated by ALK and on the subsequent detection of the phosphorylated product thus obtained. In a preferred embodiment of the invention the assay is an ELISA wherein the phosphorylated product is  
5 detected by immunochemical reactions.

In order to generate a selective ALK substrate, peptides reproducing the sequence of ALK activation loop (aa 1274-1294: ALK\_HUMAN, Q9UM73, swiss-PROT) were synthesized and tested. The peptides  
ARDIYRASFFRKGGCAMLVK (SEQ ID N. 1) and  
10 ARDIYRASYYRKGGCAMLVK (SEQ ID N.2) were particularly effective as ALK substrates showing a phosphorylation degree higher than that of polyGlu/Tyr, a random polymer which is known to be a good substrate for most tyrosine kinases. The first object of the invention is therefore a peptide having an amino acid sequence selected from SEQ ID N. 1 and SEQ ID N. 2.  
15 A further object of the invention is a method for detecting ALK tyrosine kinase activity, which essentially comprises the steps of:

- i) incubating the ALK protein or a functional derivative thereof with a peptide selected from SEQ ID N. 1 and 2 in conditions suitable for phosphorylation of the peptide;
- 20 ii) detecting the phosphorylated peptide thus formed.

As used herein, "ALK functional derivative" means any modified form of ALK protein, for example a truncated or conjugated form or a fragment thereof, which maintains the catalytic activity of unmodified ALK. The functional derivative should preferably contain the entire catalytic domain of  
25 ALK spanning residues 1116-1392 of ALK sequence (Q9UM73). The portion of ALK protein stretching from residue Leu<sup>1073</sup> to Ala<sup>1459</sup> is preferably used. When produced by recombinant gene technology using the baculovirus-based expression system, this ALK fragment shows a correct folding (confirmed by

CD spectra) and an effective catalytic activity.

In addition, a preparation containing a constitutively active form of ALK may be used instead of the purified protein or its functional derivative. Such a preparation is preferably a cell lysate, which could be used to a) assess  
5 for ALK activity in whole cells expressing ALK or ALK fusion proteins and b) assess for the effect of potential inhibitors on ALK within a cell, by treating ALK-expressing cells with different compounds and then assaying cell lysate for ALK kinase activity.

To detect the phosphorylated peptide in step ii), the phosphorylation  
10 reaction of step i) can be performed in the presence of a radioactive reagent, such as  $[\gamma^{32}\text{P}]\text{ATP}$  or  $[\gamma^{33}\text{P}]\text{ATP}$ , whereby a radioactive product is formed which is easily detected by radiometric techniques. Alternatively, an immunoreaction can be carried out, which comprises the formation of an immune complex between the phosphorylated peptide and an anti-  
15 phosphotyrosine antibody. This anti-phosphotyrosine antibody can be radioactively labelled or conjugated to a reporter enzyme, such as horse-radish peroxidase, beta-galactosidase, alkaline phosphatase or glucose oxidase, thereby allowing the direct detection of the antibody and thus of the phosphorylated peptide by measuring the specific radioactivity or enzymatic  
20 activity. As a further alternative, the anti-phosphotyrosine antibody may be detected indirectly by a second antibody recognizing the anti-phosphotyrosine antibody and carrying a radioactive label or a reporter enzyme, or by an enzyme conjugated streptavidin which recognises a biotin label on the anti-phosphotyrosine antibody.

25 In a preferred embodiment the invention provides a method for detecting ALK tyrosine kinase activity which comprises the steps of:

- a) adhering a peptide of SEQ ID N: 1 or 2 to a solid phase;
- b) incubating the solid phase with an ALK fragment extending from

Leu<sup>1073</sup> to Ala<sup>1459</sup> in suitable conditions for tyrosine phosphorylation;

- c) washing the solid phase;
- d) incubating the solid phase with an anti-phosphotyrosine antibody (primary antibody) in conditions suitable for antigen-antibody binding;
- e) washing the solid phase;
- f) incubating the solid phase with an enzyme-conjugated antibody (secondary antibody) recognizing the anti-phosphotyrosine primary antibody in conditions suitable for binding of primary antibody to secondary antibody, so that a ternary immune complex is formed;
- g) washing the solid phase;
- h) measuring the enzymatic activity of the immune complex wherein the measured activity is proportional to the amount of tyrosine-phosphorylation.

This assay, which is properly an ELISA-based kinase assay, utilizes enzyme-antibody conjugates. The conjugated enzyme cleaves a substrate to generate a colored reaction product that can be detected spectrophotometrically by measuring the absorbance of the colored solution, which is proportional to the amount of phosphotyrosines.

Solid phases or supports which can be used according to the invention comprise plastic material (reaction plates, wells, vials), polystyrene, latex and magnetic beads, colloidal metal particles, glass and/or silica surfaces and others.

The ELISA-based ALK assay is preferably used for the screening of compounds that modulate the tyrosine-phosphorylating activity of ALK, in particular for screening ALK inhibitors.


In a preferred embodiment the invention is therefore directed to a

method for the identification of compounds that modulate ALK tyrosine-kinase activity, wherein the ALK assay described above is carried out in the presence of a candidate compound or of a compound known to stimulate or inhibit ALK tyrosine kinase activity (control). Specifically, the method for identifying compounds that modulate ALK tyrosine-kinase activity comprises the steps of

- i) incubating ALK protein or a functional derivative thereof with a peptide selected from SEQ ID N. 1 or 2, preferably SEQ ID N. 1, in the presence of a candidate compound, in conditions suitable for phosphorylation of the peptide;
- ii) quantifying the phosphorylated peptide thus formed;

Compounds that inhibit ALK activity are selected as potential therapeutic agents for use in the treatment of ALK-related tumors, such as anaplastic large cell lymphomas and non-Hodgkin lymphomas. The ALK-modulating activity of a candidate compound can be compared to that of a reference compound (control), which is assayed under the same conditions as the candidate compound. Staurosporine (Meggio F. et al, Eyr. J. Biochem. (1995) 234, 317-322), which proved effective as an ALK inhibitor at 0.2 – 0.3  $\mu$ M concentration, can be used as a positive control when screening other compounds for ALK inhibitory activity.

Studies of molecular modeling gave important indications as to the pattern of substitutions required to improve staurosporine affinity and specificity towards ALK protein. The most effective structures are represented in the following general formula (I):

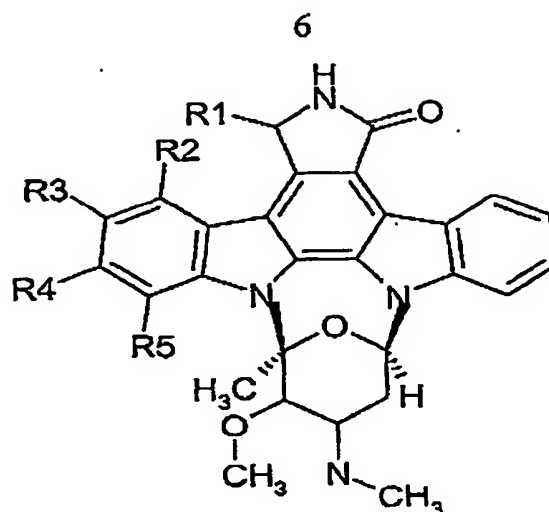

Wherein R1 and R2, independently of one another, are selected from halogen, preferably chlorine, phenyl or C1-C3 alkyl optionally substituted with one or more halogens; R3 is hydroxyl; R4 is hydroxyl or hydroxymethyl; R5 is C1-C3 alkyl, optionally halo-substituted, or benzyl.

- 5 The compounds of formula (I) antagonize the binding of ATP to the tyrosine kinase domain of ALK within oncogenic fusion proteins such as NPM-ALK, ATIC-ALK, CTLIC-ALK, TFG-ALK or other sporadic fusions with tropomyosins. Accordingly, they can be used in the preparation of a medicament for the treatment of ALK-related tumors, especially anaplastic
- 10 large cell lymphomas and non-Hodgkin lymphomas.

In a different embodiment of the invention, compounds (I) are used in an ALK-assay for the screening of ALK-inhibitors, as herein disclosed. High-throughput screenings can also be implemented using the ALK assay of the invention.

- 15 According to a further aspect, the invention is directed to a kit for carrying out the ALK assay described above. The kit consists of a packaged combination of reagents in predetermined ratios. Typically the kit will contain a peptide of SEQ ID N: 1 or 2, optionally immobilized on a solid phase, the anti-phosphotyrosine antibody and, if necessary, an antibody labelled with an
- 20 enzymatic or radioactive label. In addition the kit may contain reagents for



colorimetric or radiometric reactions, buffers, controls such as staurosporine or a derivative thereof, diluents, detergents, stabilizers and any other component useful for setting up and carrying out the assay. The relative amounts of the various reagents may be varied to optimize the sensitivity of the assay. Particularly the reagents may be provided in solid or liquid preparations, such as solution, suspension, dispersion, dry powders, lyophilized preparation. The kit components are supplied in single or separate containers. The kit may also include instructions for carrying out the assay.

### DESCRIPTION OF THE FIGURES

10 **Figure 1: purification and activity of recombinant His-tagged ALK**

A) Silver stained 10% SDS-PAGE gel showing marker (M), dialysed and pooled DEAE column fractions loaded onto the HiTrap column (load), HiTrap column fractions (numbers indicate fraction number). B) Radioactive autophosphorylation assay of purified rALK. Lane 1: autoradiograph of  $^{32}\text{P}$ -labelled rALK. Lane 2: silver staining of the same sample as in Lane 1.

**Figure 2: kinetics for ALK kinase with peptides**

**Figure 3: Detection of rALK activity using the ELISA based kinase assay**

A) Phosphorylation of peptide SEQ ID N.2 by purified rALK. An ELISA kinase reaction was performed with or without 0.5  $\mu\text{g}$  purified rALK, 15  $\mu\text{g}$  peptide SEQ. ID N.2 or no peptide, at 30°C for 30 mins. The graph shows absorbance at 450 nm. B) The kinetics of ALK peptide substrate phosphorylation. An ELISA kinase reaction was performed using 206  $\mu\text{M}$  peptide SEQ. ID N.2 or 42  $\mu\text{M}$  peptide SEQ. ID N.1 with 0.1  $\mu\text{g}$  of purified rALK. The reaction was stopped by adding EDTA after 0, 0.5, 2, 5, 10 and 15 mins.

**Figure 4: The effect of peptide substrate concentration on level of phosphorylation**

The ELISA assay was performed in the presence of 0.2  $\mu$ g purified rALK and 0 - 105  $\mu$ M peptide SEQ ID N.1 (A) or 0 - 315  $\mu$ M peptide SEQ ID N.2 (B), at 30°C for 15 mins.

5 **Figure 5: The effect of rALK concentration on substrate phosphorylation**

An ELISA assay was performed using 206  $\mu$ M peptide SEQ ID N.2 or 42  $\mu$ M peptide SEQ ID N.1 and rALK ranging from 0 – 225 ng, at 30°C for 15 mins.

**Figure 6: Inhibition of rALK activity by Staurosporine**

10 An ELISA assay was performed using 42  $\mu$ M SEQ ID N.1, 20 ng rALK, at 30°C for 15 mins in the presence of 0 – 5  $\mu$ M staurosporine or an equivalent volume of solvent, DMSO. Graph shows A450 normalised to the untreated control.

## EXPERIMENTAL RESULTS

### 15 Peptide synthesis

Peptides were synthesised by automatic solid phase synthesis utilising 9-fluorenylmethoxycarbonyl (Fmoc) chemistry with an Applied Biosystems Model 431A synthesiser on 4-hydroxymethyl-copolystyrene-1% divinylbenzene-resin (0.95 mmol/g, 0.05 mmol).

20 Fmoc amino acids (0.5 mmol) were activated by 2-(1H-benzotriazol-1-yl)-1,1,3,3-tetramethyluronium hexafluorophosphate (HBTU) (1 eq) and 1-hydroxybenzotriazole (HOBt) (1 eq) in the presence of DIPEA (2 eq).

25 Peptidyl-resins (100 mg) were cleaved and deprotected in a mixture containing 5 ml of trifluoroacetic acid, 0.375 gr of phenol, 0.125  $\mu$ l of 1,2-ethanedithiol, 0.250  $\mu$ l of thioanisole, and 0.25  $\mu$ l of water for 2 hours. The reaction mixture was filtered onto a tube containing ethyl ether cooled at 0°C. Precipitated peptides were separated by centrifugation and washed with fresh ether.

Crude peptides (50-100 mg in 10 ml of water) were pumped onto a preparative RP-column (prepNova-Pak HR C18, 6  $\mu$ m, 25x10 mm, Waters) and eluted with a linear gradient of 10% to 45% acetonitrile at 12 ml/min.

The purity of peptides was > 90% as judged by analytical RP-HPLC on a 5  $\mu$ m, C18 Symmetry300 column, 4.6x250 mm (Waters) using a linear gradient of acetonitrile in 0.1% trifluoroacetic acid at 1 ml/min. Peptide molecular weights were confirmed by mass spectroscopy using a matrix-assisted laser-desorption ionisation time-of-flight (MALDI-Tof) spectrometer (Maldi-1; Kratos-Schimidzu, Manchester, UK).

#### 10      **Production and purification of recombinant ALK**

A portion of ALK, stretching from amino acids Leu<sup>1073</sup> to Ala<sup>1459</sup>, has been cloned into the baculovirus transfer vector, pBlueBacHis2C (Invitrogen). Using this vector and the MaxBac® 2.0 baculovirus expression system (Invitrogen) we have produced recombinant baculovirus and expressed recombinant ALK (rALK) protein in Sf9 (*Spodoptera frugiperda*) insect cells. This protein has a theoretical molecular weight of 49 kDa and contains the predicted catalytic domain of ALK (Ile<sup>1116</sup> to Val<sup>1392</sup>) fused to a 6-histidine-tag (His-tag). The autophosphorylation activity of His-tagged rALK, which is indicative of correct folding, has been verified by anti-phosphotyrosine immunoblotting of whole cell lysates and an *in vitro* radioactive kinase assay.

To produce protein for purification, Sf9 cells were infected with recombinant virus at a multiplicity of infection (MOI) of 5. Cultures were incubated for 3 days at 27°C and then harvested by centrifuging at 400 g for 10 minutes at 4°C. Cell pellets were stored at -80°C until use. To lyse the cells, cell pellets were resuspended in hypotonic buffer (Buffer A) containing 50 mM Tris-HCl pH 8, 20 mM NaCl and protease inhibitors (Pepstatin, Benzamidine, Leupeptin and Aprotinin). After 30 minutes incubation on ice, the cell suspension was centrifuged at 1500 g for 15 minutes at 4°C and the

supernatant was filtered through a 0.45  $\mu\text{m}$  filter. The filtrate was then loaded onto an 80 ml anion exchange Fast Flow DEAE-sepharose (Sigma) column in the presence of Buffer A at a flow of 1.5 ml/min using the AKTA FPLC system (Amersham-Pharmacia Biotech). Bound proteins were eluted in gradient of NaCl ranging from 20 to 200 mM. Fractions were analysed for the presence of ALK protein by immunoblotting. Positive fractions were pooled and dialysed for 3 h at 4°C in 1 litre of native binding buffer (20 mM sodium phosphate pH 7.8, 500 mM NaCl and PI), changing the buffer every hour. Dialysed fractions were then loaded onto a HiTrap™-nickel affinity chromatography column (Amersham-Pharmacia Biotech), previously equilibrated with native binding buffer supplemented with 50 mM imidazole and 20 mM  $\beta$ -mercaptoethanol. A flow of 0.5 - 1 ml/min was applied to the column. After washing the column with native binding buffer bound proteins were eluted using a linear gradient of imidazole ranging from 50 to 200 mM. Fractions were then analysed for the presence and purity of rALK by SDS-PAGE and silver staining. Figure 1A shows an example purification of rALK.

Finally, positive fractions were dialysed in 50 mM Tris pH 7.4, 100 mM NaCl, 10% glycerol, 20 mM  $\beta$ -mercaptoethanol and protease inhibitors, and stored at -80°C until use. To verify that purified rALK was active we have performed a radioactive autophosphorylation assay. A reaction mix containing 25  $\mu\text{l}$  of a rALK positive fraction, 6  $\mu\text{M}$  cold ATP, 1.2 mM DTT, 10  $\mu\text{Ci}$  [ $\gamma$ - $^{32}\text{P}$ ]ATP, 25 mM Hepes pH 7.5, 10 mM  $\text{MgCl}_2$  and 10 mM  $\text{MnCl}_2$  was incubated at 30°C for 30 minutes. The reaction was stopped by adding Laemmli buffer and heating at 95°C for 5 minutes. The samples were resolved on a 10% SDS-PAGE gel and transferred to an Immobilon™-P membrane. Radioactively labelled proteins were visualized by autoradiography as shown in Figure 1B.

### Peptide phosphorylation (radioactive assay)

Peptides and the random polymer polyGlu/Tyr (1:4) were phosphorylated in 30  $\mu$ l of a medium containing 50 mM Tris/HCl, pH 7.5, 5 mM  $\text{MnCl}_2$ , 10  $\mu$ M Na-vanadate, 30  $\mu$ M  $[\gamma^{32}\text{P}]\text{ATP}$  or  $[\gamma^{33}\text{P}]\text{ATP}$  (specific activity 1000 cpm/pmol) and 10 units of rALK [one unit was defined as the amount of enzyme transferring 1 pmol phosphate per min to polyGluTyr (0.1 mg/ml)]. The reactions were terminated after 10 min of incubation at 30°C by spotting 25  $\mu$ l of the mixture onto P81 phosphocellulose papers, which were processed as described elsewhere (1). Kinetic constants were determined by GraphPad Prism software fitting the data directly to the Michaelis-Menten equation using nonlinear regression.

### Results

The peptide derived from ALK reference sequence, containing Tyr-1278 was a good substrate of rALK, whereas peptides containing either Tyr-1282 or Tyr-1283 were slightly affected by the enzyme (see Fig. 2). The Tyr-1278-derivative was phosphorylated with an efficiency even higher than that displayed by the peptide bearing three Tyr residues (Table I). To confirm that Tyr-1278-derivative is an optimal peptide substrate for ALK catalytic domain, we compared its phosphorylation degree with that displayed by polyGlu/Tyr (1:4), a random polymer that is a very good substrate for most tyrosine kinases. Table II shows that Tyr-1278-derivative is a better ALK substrate than polyGlu/Tyr.

TABLE I

Kinetic constants for rALK with synthetic peptides.

PEPTIDE	$V_{\max}$ (pmol/min)	$K_m$ ( $\mu$ M)	Efficiency ( $V_{\max}/K_m$ )
ARDIYRASYYRKGGCAML PVK	99.5	90.5	1.1
ARDIYRASFFRKGGCAML PVK	186.3	109.4	1.7

TABLE II

Phosphorylation rates of model substrates by ALK catalytic domain.

Poly(Glu/Tyr) and peptide concentrations were 0.1 mg/ml and 400  $\mu$ M, respectively. Enzyme concentration was 10 units. Reported values represent the means for three separate experiments. S.E.M. values were always less than 14%.

SUBSTRATE	Phosphorylation degree (pmol/min)
Poly(Glu/Tyr)	10.0
ARDIYRASFFRKGGCAMLVK	30.3

#### ELISA-based *in vitro* kinase assay

A Nunc Immuno 96 well plate was incubated overnight at 37°C with coating solution (125  $\mu$ l/ well) containing ALK peptide substrate (SEQ ID N. 1 or SEQ ID N. 2) at various concentrations in PBS. The wells were then washed with 200  $\mu$ l of wash buffer (PBS-Tween 0.05 %) and left to dry for at least 2 hours at 37°C. The kinase reaction was performed by incubating kinase buffer (25 mM Hepes pH 7.5, 5 mM  $MnCl_2$ , 5 mM  $MgCl_2$ ), 0.3 mM ATP and purified rALK at various concentrations in a total volume of 100  $\mu$ l/well at 30°C for 15 minutes. The reaction mix was then removed and wells were washed 5 times with 200  $\mu$ l of wash buffer. Phosphorylated peptide was detected using 100  $\mu$ l/well of a mouse monoclonal anti-phosphotyrosine antibody (clone 4G10 UpstateBiotech Ltd) diluted 1:500 in PBS + 4% BSA. After 30 minutes incubation at room temperature the antibody was removed and wells were washed as described above. 100  $\mu$ l of a secondary antibody (anti-mouse IgG, Horseradish Peroxidase linked whole antibody, from sheep, Amersham Pharmacia Biotech) diluted 1:1000 in PBS + 4% BSA was added to each well and the plate was incubated again for 30 minutes at room temperature before washing as above. The plate was developed using 100

μl/well TMB Substrate Solution (Endogen) and the reaction was stopped by adding an equal volume of H<sub>2</sub>SO<sub>4</sub> 0.18 M. Finally, the absorbance was read at 450 or 490 nm using an Ultrospec® 300 spectrophotometer (Amersham-Pharmacia Biotech).

## 5 *Results*

We have shown that the ELISA-based kinase assay can be used to detect phosphorylated ALK peptide substrates (SEQ ID N.1 and N.2) and hence the activity of purified rALK. In initial experiments we performed the assay in the presence and absence of both purified rALK and the peptide  
10 substrate SEQ ID N.2, in order to determine if a specific increase in absorbance at 450 nm (A<sub>450</sub>) could be observed (Figure 3A). Indeed in the presence of both rALK and substrate a dramatic increase in A<sub>450</sub> was observed. In contrast, in the absence of either peptide or rALK or both, A<sub>450</sub> was low, indicating a low level of background absorbance. To determine the  
15 kinetics of substrate phosphorylation we measured A<sub>450</sub> after various incubation times (Figure 3B). For both peptides the kinetics of the reaction was similar with maximum phosphorylation achieved after 10 minutes.

To compare the efficiency of the two peptides in the ELISA kinase assay and to determine the optimal concentration of peptide to use we  
20 performed standard curves for both peptides. Results show that increasing the amount of peptide substrate increased A<sub>450</sub> and hence rALK kinase activity. For the peptide substrate SEQ ID N.2, maximum A<sub>450</sub> and hence phosphorylation was observed at approximately 200 μM, compared with only 25 μM for peptide substrate SEQ ID N.1 (Figures 4A and B). This indicates  
25 that peptide substrate SEQ ID N.1 is more efficient as a substrate for rALK.

To observe inhibition of kinase activity it is essential that we use an appropriate concentration of enzyme, i.e. in the linear range of a standard

curve. An ELISA was performed using 206  $\mu\text{M}$  SEQ ID N.2 or 42  $\mu\text{M}$  SEQ ID N.1 and rALK ranging from 0 – 220 ng. The linear range of the curve for both substrates was approximately between 0 – 15 ng of rALK, after which the curves reached plateau (Figure 5).

- 5           The effect of the inhibitor staurosporine on rALK activity was assessed using 42  $\mu\text{M}$  peptide substrate SEQ ID N.2. A solvent control containing DMSO at corresponding concentrations was also performed. Staurosporine markedly inhibited rALK activity reaching maximum inhibition at 1  $\mu\text{M}$  and with an estimated  $\text{IC}_{50} = 300 \text{ nM}$  (Figure 6). The solvent (DMSO) alone had
- 10 no substantial effect on rALK activity.



**REFERENCES**

- [1] Glass, D.B., Masaracchia, R.A., Feramisco, J.R. and Kemp, D.E. (1978) Anal. Biochem. 87, 566-575.
- 5 [2] Till, J.H., Ablooglu, A.J., Frankel, M., Bishop, S.M., Kohanski, R.A. and Hubbard S.R. (2001) J. Biol. Chem. 276, 10049-10055.

CLAIMS

1. A method for detecting ALK tyrosine kinase activity, which comprises the following steps:

5           i) incubating the ALK protein or a functional derivative thereof with a peptide substrate selected from SEQ ID N. 1 or 2 in conditions suitable for phosphorylation of the peptide;

          ii) detecting the phosphorylated peptide.

10       2. A method according to claim 1, wherein the peptide has sequence SEQ ID N. 1.

3. A method according to claim 1, wherein purified ALK protein or an ALK-containing preparation is used.

4. A method according to claim 3, wherein said preparation is a cell lysate.

15       5. A method according to claim 1, wherein said functional derivative contains the entire catalytic domain of ALK spanning residues 1116-1392 of ALK sequence.

6. A method according to claim 5, wherein said functional derivative is a fragment of ALK protein extending from residue Leu<sup>1073</sup> to Ala<sup>1459</sup>.

20       7. A method according to claim 6, which comprises the steps of:

a) adhering a peptide of SEQ ID N. 1 or 2 to a solid phase;

b) incubating the solid phase with said ALK fragment in conditions suitable for tyrosine phosphorylation;

c) washing the solid phase;

25       d) incubating the solid phase with an anti-phosphotyrosine antibody (primary antibody) in conditions suitable for antigen-antibody binding;

e) washing the solid phase;

f) incubating the solid phase with an enzyme-conjugated antibody

(secondary antibody) recognizing the primary antibody in conditions suitable for the binding of primary and secondary antibodies, so that a ternary immune complex is formed;

g) washing the solid phase;

5 h) measuring the enzymatic activity of the immune complex wherein the measured activity is proportional to the amount of tyrosine-phosphorylation.

8. A method according to claim 7, wherein the enzyme conjugated to the antibody is Horse-Radish peroxidase.

10 9. A method according to claim 7, wherein the enzymatic activity is detected by colorimetric reaction.

10. A method according to any previous claims, for the identification of compounds that modulate ALK tyrosine-kinase activity.

11. A method according to claim 10, which comprises the steps of


15 i) incubating ALK protein or a functional derivative thereof with a peptide selected from SEQ ID N. 1 or 2 in the presence of a candidate compound (a) in conditions suitable for peptide phosphorylation;

ii) detecting the phosphorylated peptide thus formed;

20 12. A method according to claims 10-11, wherein the ALK-modulating activity of the candidate compound is compared to that of a reference compound which is assayed under the same conditions as the candidate compound.

13. A method according to claim 12, wherein the reference compound is  
25 staurosporine.

14. A method according to claim 12, wherein the reference compound is a staurosporine derivative of general formula (I):

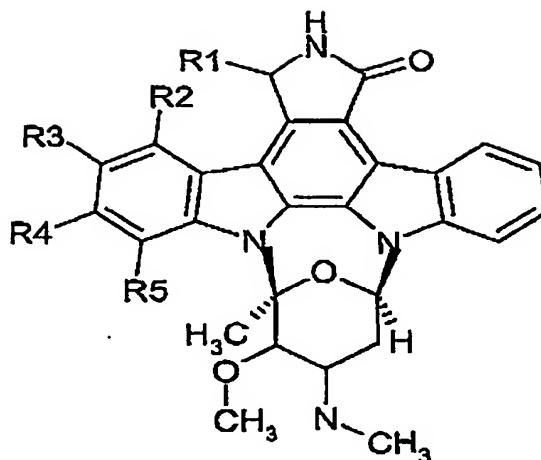

Wherein R1 and R2, independently of one another, are selected from halogen, preferably chlorine, phenyl or C1-C3 alkyl optionally substituted with one or more halogens; R3 is hydroxyl; R4 is hydroxyl or hydroxymethyl; R5 is C1-C3 alkyl, optionally halo-substituted, or benzyl.

- 5 15. A peptide useful as ALK substrate selected from SEQ ID N. 1 or 2.
16. A peptide according to claim 15, which is SEQ ID N. 1.
17. The use of a peptide according to claim 15 or 16 for the determination of ALK tyrosine-kinase activity.
18. The use of a compound of formula (I), as per claim 14, for the preparation of a medicament for the treatment of ALK-related tumors, especially anaplastic large cell lymphomas and non-Hodgkin lymphomas.
- 10 19. A kit for detecting ALK tyrosine-kinase activity according to claims 1-14, which comprises a peptide of SEQ ID N: 1 or 2 and an anti-phosphotyrosine antibody.
- 15 20. A kit according to claim 19, containing an additional component selected from reagents for colorimetric reactions, buffers, diluents, detergents, stabilizers, staurosporine or a derivative thereof as per claim 14.

## SEQUENCE LISTING

<110> ISTITUTO NAZIONALE PER LO STUDIO E LA CURA DEI TUMORI

<120> ANAPLASTIC LYMPHOMA KINASE ASSAY, REAGENTS AND COMPOSITIONS THEREOF

<130> 1206EUR

<160> 2

<170> PatentIn version 3.1

<210> 1

<211> 21

<212> PRT

<213> synthetic oligonucleotide

<400> 1

Ala Arg Asp Ile Tyr Arg Ala Ser Phe Phe Arg Lys Gly Gly Cys Ala  
1 5 10 15

Met Leu Pro Val Lys  
20

<210> 2

<211> 21

<212> PRT

<213> synthetic oligonucleotide

<400> 2


Ala Arg Asp Ile Tyr Arg Ala Ser Tyr Tyr Arg Lys Gly Gly Cys Ala  
1 5 10 15

Met Leu Pro Val Lys  
20


Sheet 1/6

Figure 1

A.

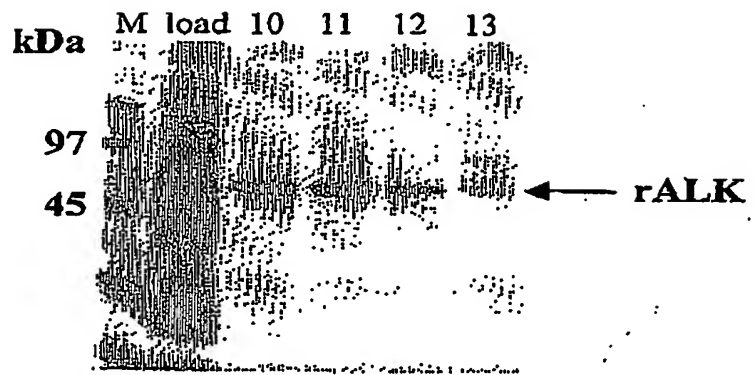


B.

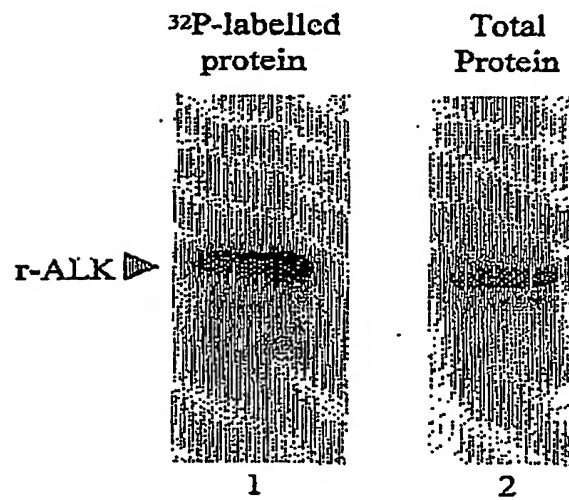


## Sheet 2/6

Figure 2

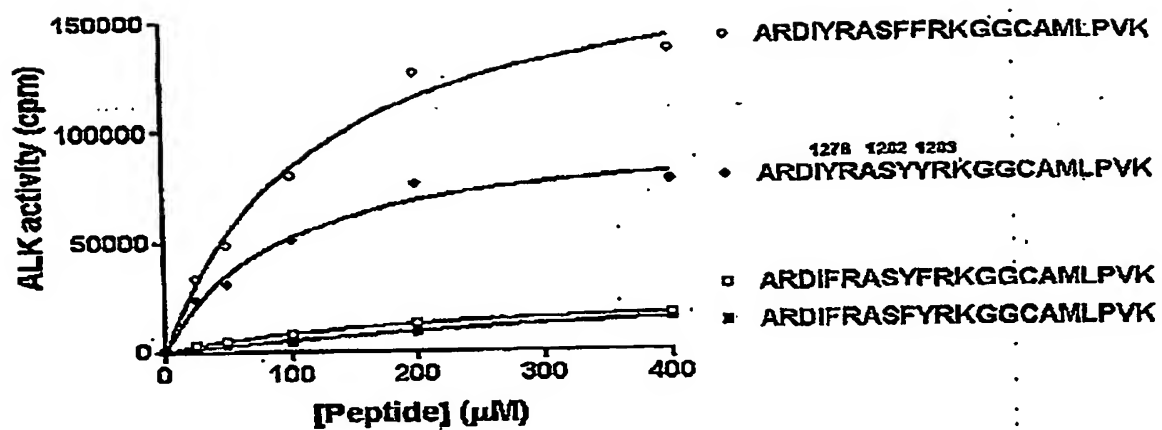

## Sheet 3/6

### Figure 3

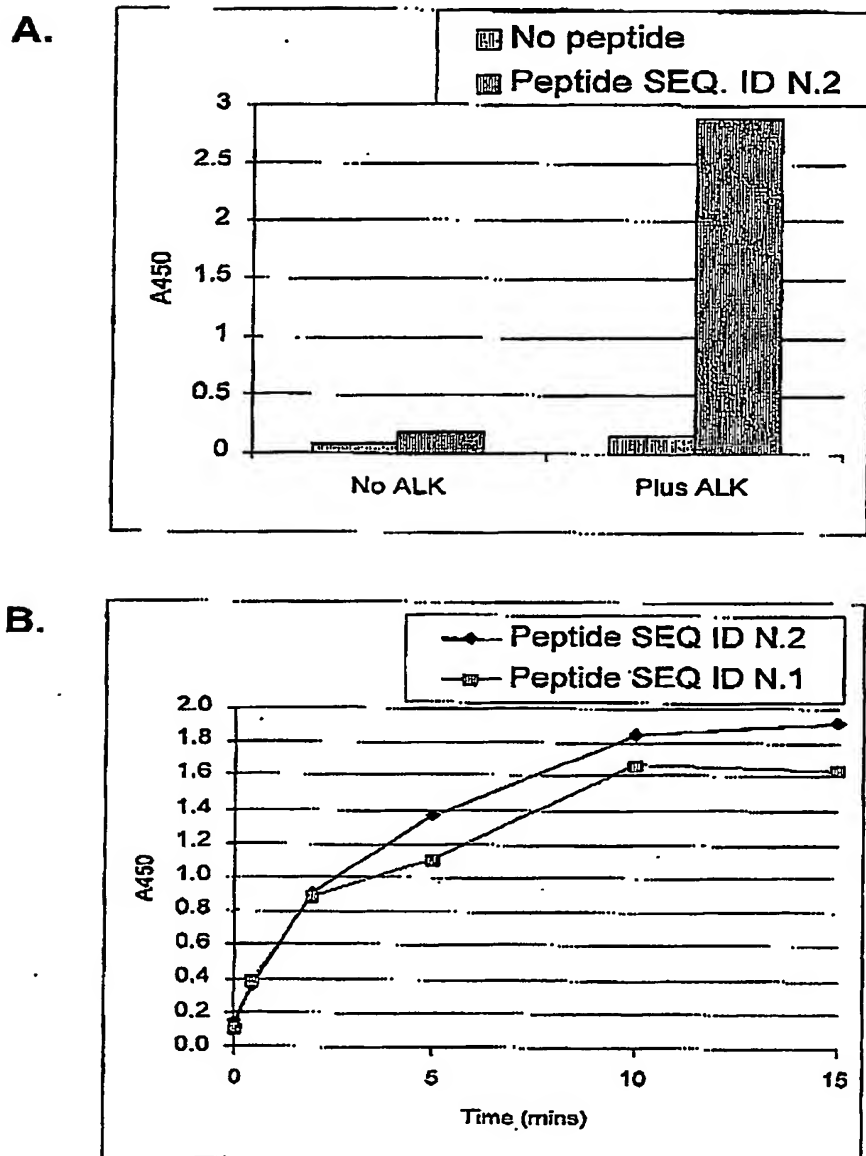




# Sheet 4/6

## Figure 4

A.

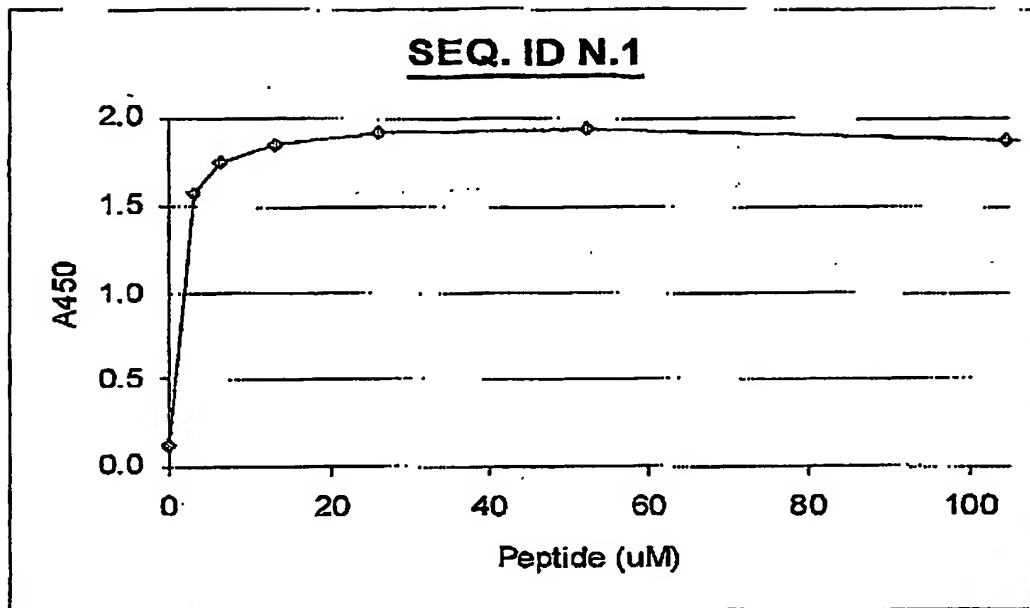

B.

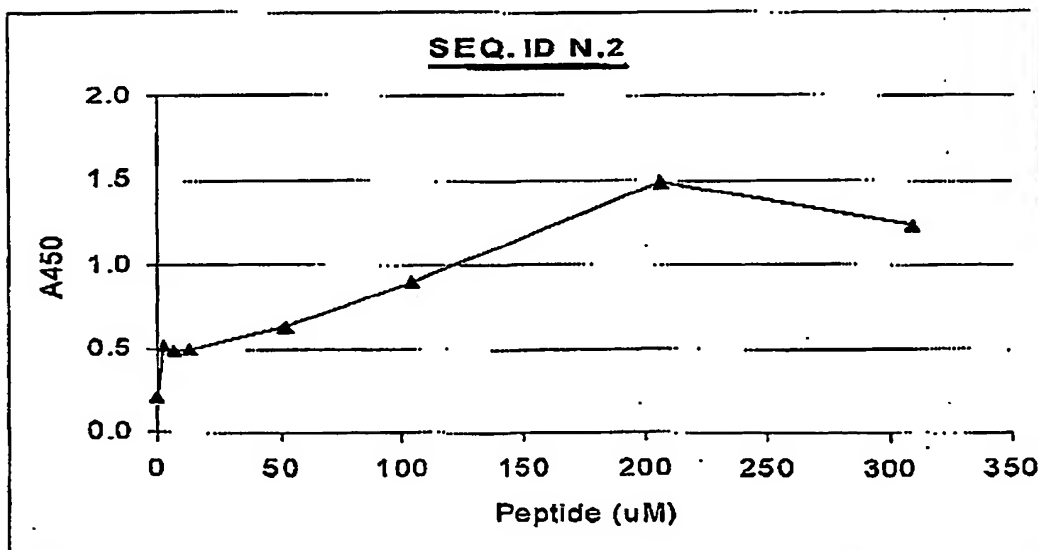


Figure 5


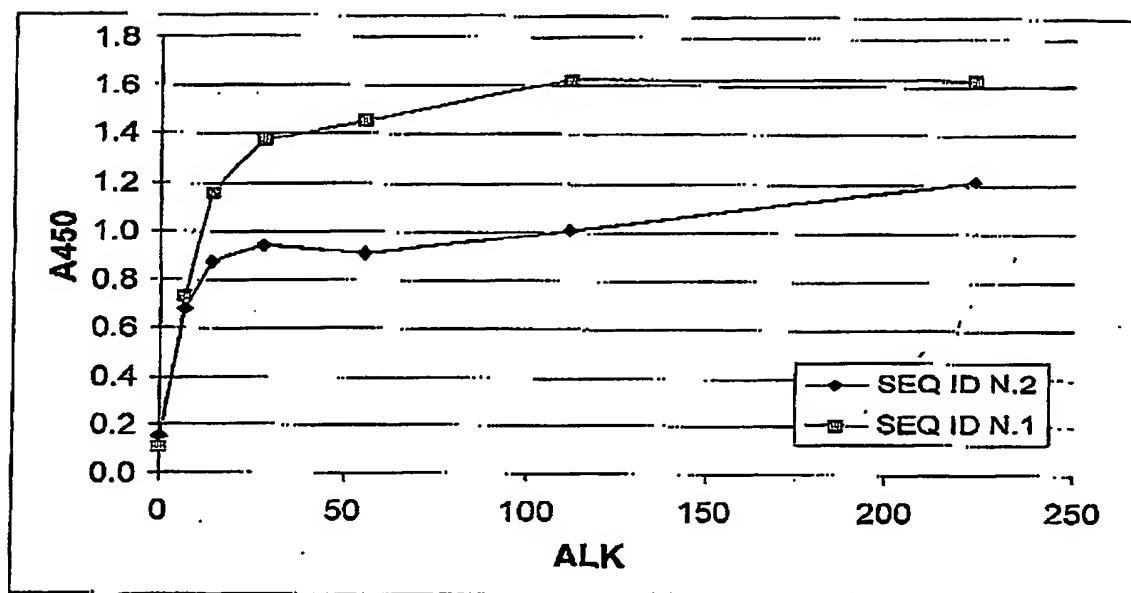


Figure 6

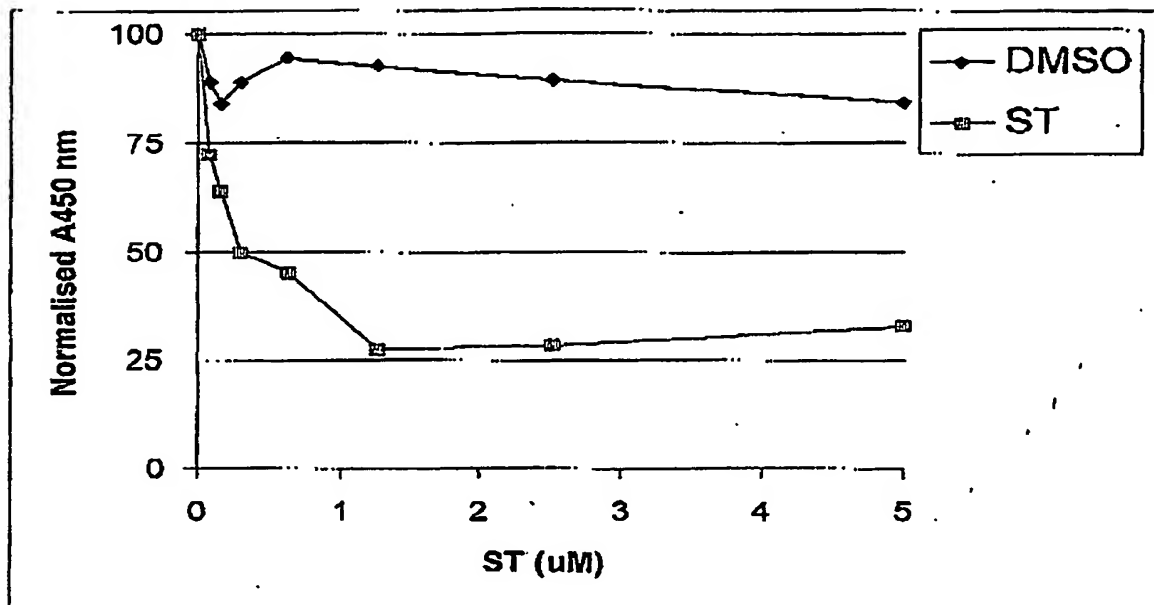

**ABSTRACT****ANAPLASTIC LYMPHOMA KINASE ASSAY, REAGENTS AND COMPOSITIONS THEREOF**

5

Disclosed are methods for detecting ALK tyrosine kinase activity and for the identification of compounds that modulate ALK activity, peptide substrates, reagents and compositions thereof.

EP 04 2185



This Page is inserted by IFW Indexing and Scanning Operations and is not part of the Official Record

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☐ BLURED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLORED OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REPERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images problems checked, please do not report the problems to the IFW Image Problem Mailbox**